


## SofTouch<sup>™</sup> Elite Bellyboard



Shown with Prone Headrest



15 mm headrest

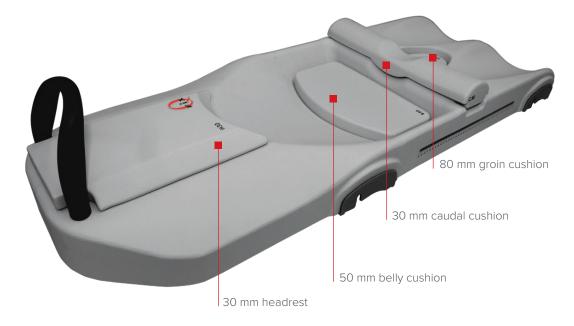

## SofTouch™ Elite Bellyboard

The SofTouch™ Elite Bellyboard is an adjustable stabilization platform designed to position the lower bowel away from the treatment field during prostate radiotherapy. The removable caudal cushion and separate belly support cushions allow for custom patient setup. The SofTouch Elite Bellyboard's unique semi-firm radiolucent foam provides repeatable positioning and is contoured to relieve pressure point stress and improve patient comfort. Multiple cushion configurations provide optimal anatomical placement of both external and internal anatomy for improved clinical outcomes.

## Features:

- Three contoured and adjustable headrests allow for improved patient comfort
- · Provides multiple indexing locations to ensure precise, repeatable positioning
- Ergonomic hand grip relieves stress on patient's wrists
- Removable cushions allow for custom patient setup
- Contoured to relieve pressure point stress
- Scales to aid in laser alignment
- Lightweight Easy to clean

RT-4464-02 SofTouch™ Elite Bellyboard



